

RBRC 00067

Strain name BALB/cA-TgH(Tnc)

Primer	primer_name	primer_seq	length	
1	TN-G1-1	GGTACCTGATTCCGAAGTGCATTGTCACG	29	Depositor's primer
2	TN-150R (~3'TN)	AAGCTAGAAAGATGCCTGGC	20	Depositor's primer
3	ZA-6	ACGGCAAACGACTGTCCTGGCCGTAACCGACC	32	Depositor's primer

Reaction Components


Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		1.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	120*	*Repeat for 30 cycles

Standard PCR

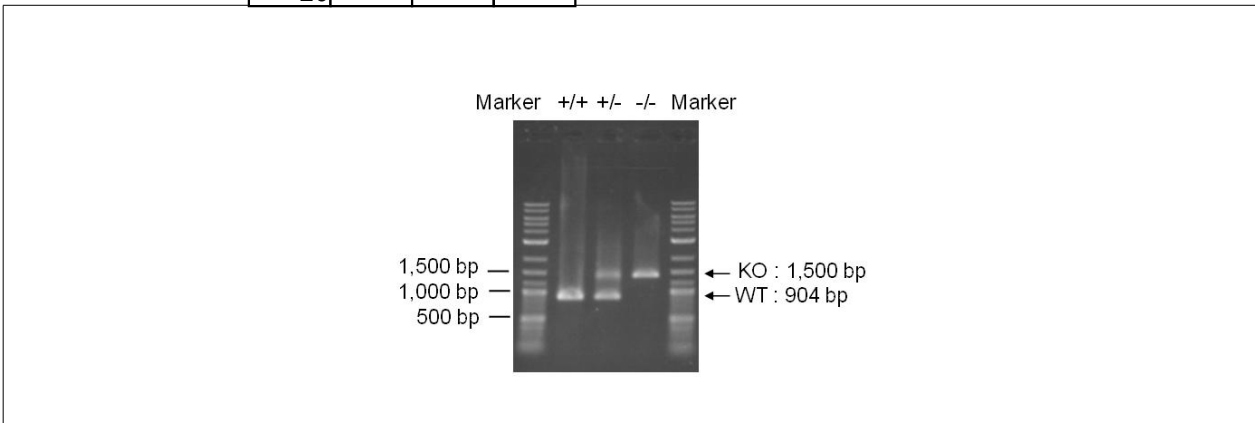

Marker
1 kb DNA ladder (Invitrogen, Cat. 15615-024)

Gel
1% agarose in 1 X TAE (100 V, 30 min)

PCR products	primer_set	product_size	product
primer 1	primer 2	904 bp	WT
primer 1	primer 3	1,500 bp	KO

Comment

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.