

RBRC 02095 Strain name B6.B6D2-Tg(Th-EGFP)21-31Koba

Primer	primer_name	primer_seq	length	
1	FW_Primer	CCTGTGACAGTGGATGCAATT	21	RIKEN BRC
2		CTTGTACAGCTCGTCCATGCCGAG	24	RIKEN BRC
3	mLC3ex3GT	TGAGCGAGCTCATCAAGATAATCAGGT	27	Mol. Biol. Cell., 15, 1101-1111 (2004).
4	mLC3ex4AG	GTTAGCATTGAGCTGCAAGCGCCGTCT	27	Mol. Biol. Cell., 15, 1101-1111 (2004).


Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	120*	*Repeat for 30 cycles

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/μl	1.0			
primer 2	10 pmol/μl	1.0			
primer 3	10 pmol/μl	1.0			
primer 4	10 pmol/μl	1.0			
DNA		1.0			
Taq	1 U/μl	0.2			
total		20			

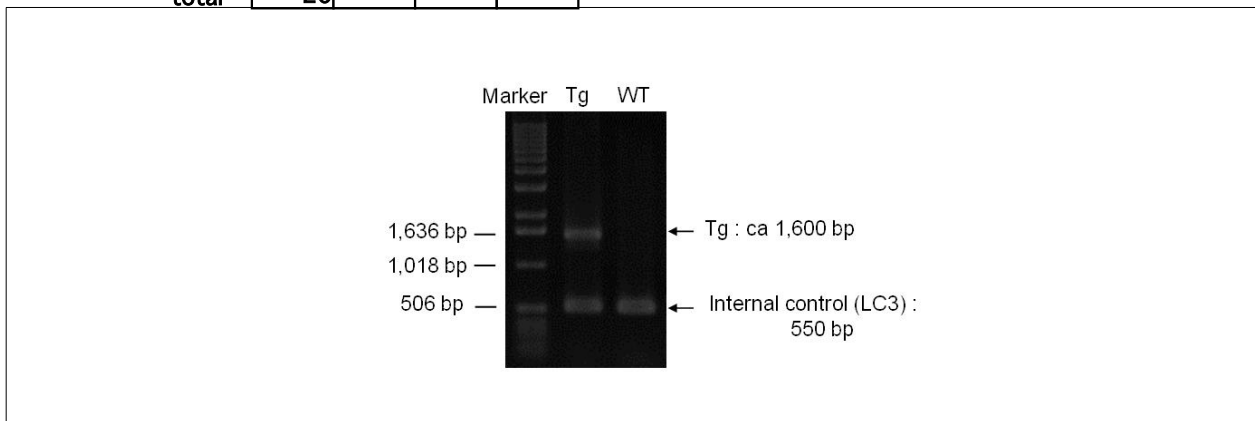

Marker
1 kb DNA ladder (Invitrogen, Cat. 15615-024)

Gel
1% agarose in 1 X TAE (100 V, 30 min)

PCR products				
	primer_set	product_size	product	
	primer 1	primer 2	ca 1600 bp	Tg
	primer 3	primer 4	550 bp	Internal control (LC3)

Comment Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.