

RBRC 02850

Strain name FVB/N- *Tg(Alb-Ppp1r15a)14Dron*

Primer	primer_name	primer_seq	length	Reference
1	promoter	ACAGCTCCAGATGGCAAACATAC	23	RIKEN BRC
2	haGADD34.6AS	CTGTGTCCTCTCCAGGGCAACAGAC	25	Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol

Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions


step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Reagent	Conc	Reaction 1	2	3	4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/μl	1.0			
primer 2	10 pmol/μl	1.0			
primer 3	10 pmol/μl	1.0			
primer 4	10 pmol/μl	1.0			
DNA		1.0			
Taq	1 U/μl	0.2			
total		20			

(μL)

Marker Gene Ladder 100 (Wako, Cat. 316-06951)

Gel 2% agarose in 1 X TAE (100 V, 30 min)

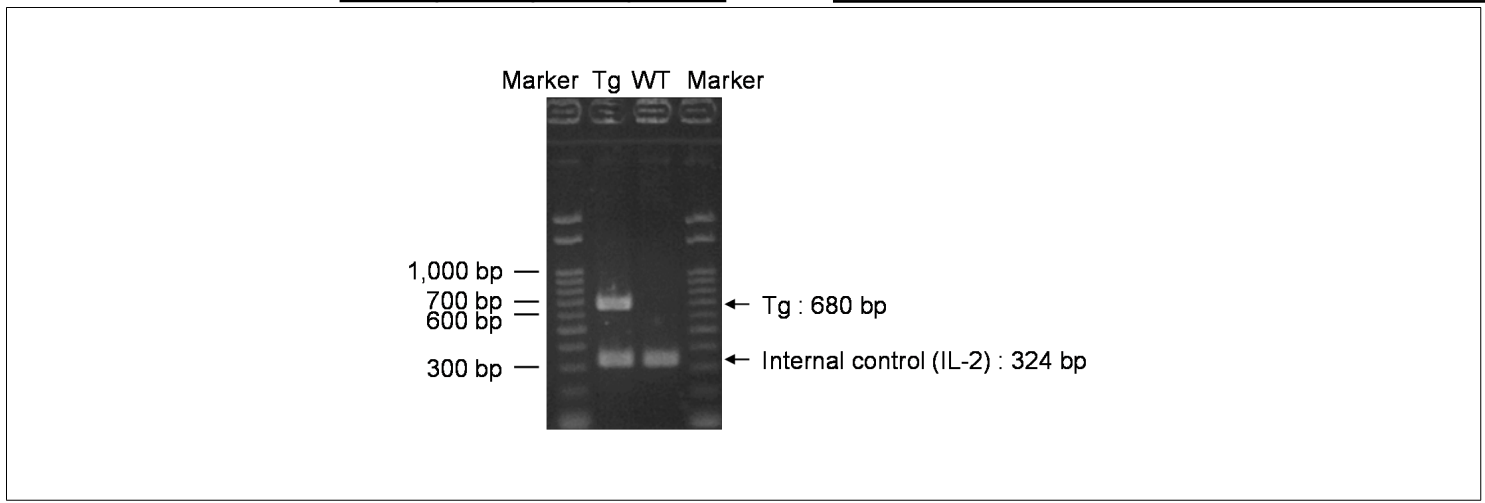

primer_set	product_size	product (allele type)
primer 1	primer 2	680 bp
primer 3	primer 4	324 bp
		Tg
		Internal control (IL-2)
		Non-Tg: 324 bp
		Tg: 324 bp, 680 bp

Comment Standard PCR

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.

