

RBRC 03513

Strain name *Isx(lacZ)* (Acc.No. CDB0475K)

Primer	primer_name	primer_seq	length	
1		AAGTGTTCCCTCCTAGCAGCA	20	J Biol Chem., 283(8):4905-4911 (2008).
2		CAATGACCTGAGCCTATGTG	20	J Biol Chem., 283(8):4905-4911 (2008).
3		ATTCAGGCTGCGCAACTGTT	20	J Biol Chem., 283(8):4905-4911 (2008).


Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		1.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

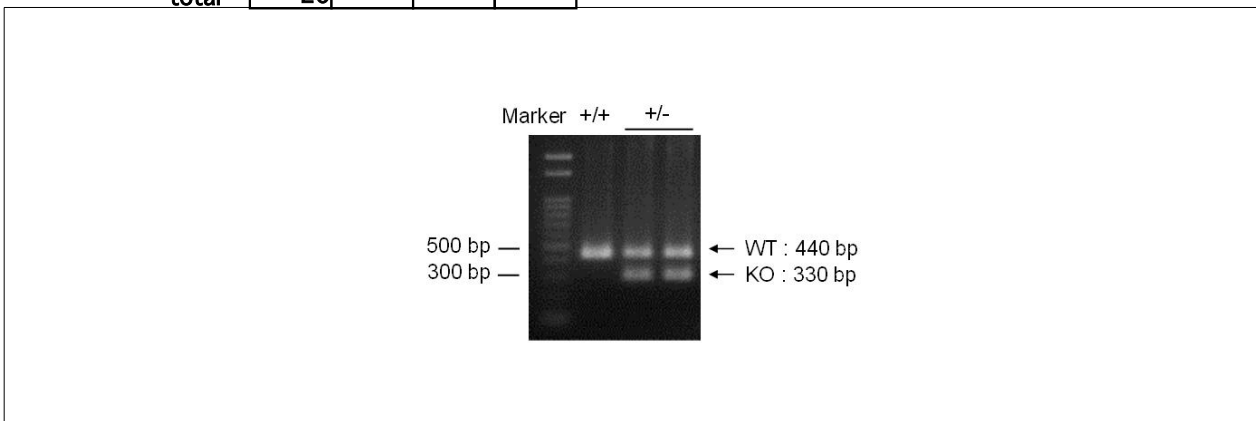

Marker
Gene Ladder 100 (Wako, Cat. 316-06951)

Gel
2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1 primer 2	440 bp	WT
primer 1 primer 3	300 bp	KO

Comment Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.