

PCR protocol

RBRC 05154

Strain name B6;129S6-Gt(ROSA)26Sor<tm1.1(CAG-EGFP/TeNT)Imayo>

Primer	primer_name	primer_seq	length	
1	R26F2	AAAGTCGCTCTGAGTTGTTAT	21	Depositor's primer
2	oIMR0316	GGAGCGGGAGAAATGGATATG	21	JAX Genotyping Protocol
3	oIMR0315	GCGAAGAGTTTGTCTCAACC	21	JAX Genotyping Protocol

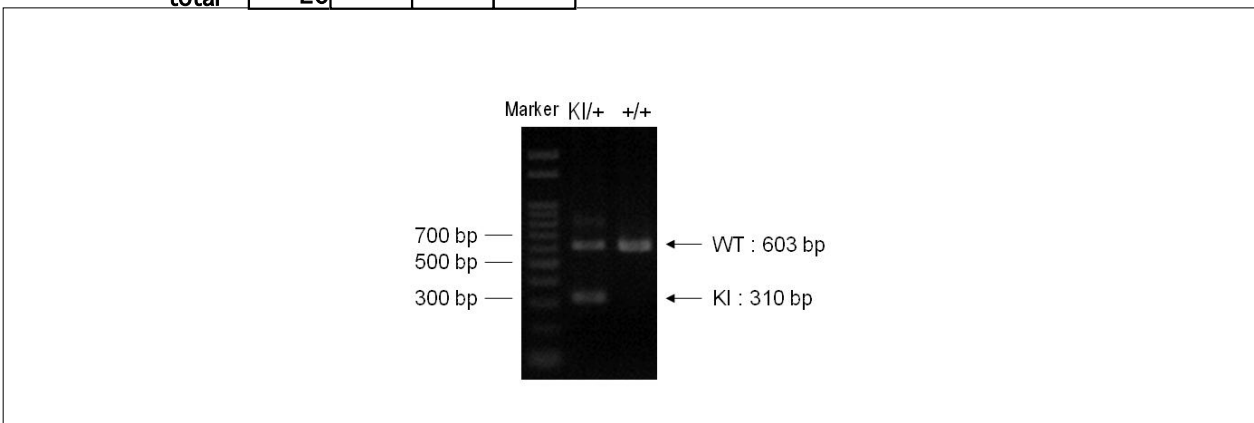
Reaction Components

Enzyme	TaKaRa Ex Taq (Cat. No. RR001A)
--------	---------------------------------

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		12.3			
PCR buffer	10 X	2.0			
dNTP	2.5 mM	1.6			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	5 U/ μ l	0.1			
total		20			

step	temp (°C)	time (sec)	note
1	94	180	
2	94	30*	
3	60	30*	
4	72	60*	*Repeat for 30 cycles
5	72	300	


Marker

Gene Ladder 100 (Wako, Cat. 316-06951)

Gel

2% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1 primer 2	603 bp	WT
primer 1 primer 3	ca 310 bp	KI

Comment