

RBRC 05250 Strain name B6;B6D2-Tg(Pou2f3-Ascl1-ires-nEGFP)2Abek

Primer	primer_name	primer_seq	length	
1		CTCCAACGACTTGA ACTCTATGGC	24	RIKEN BRC
2	EGFP-R1	CTTGTACAGCTCGTCCATGCCGAGAGTGA	29	Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol


Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		0.6			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/μl	1.0			
primer 2	10 pmol/μl	1.0			
primer 3	10 pmol/μl	1.0			
primer 4	10 pmol/μl	1.0			
DNA		1.0			
Taq	1 U/μl	0.4			
total		20			

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	90*	*Repeat for 30 cycles

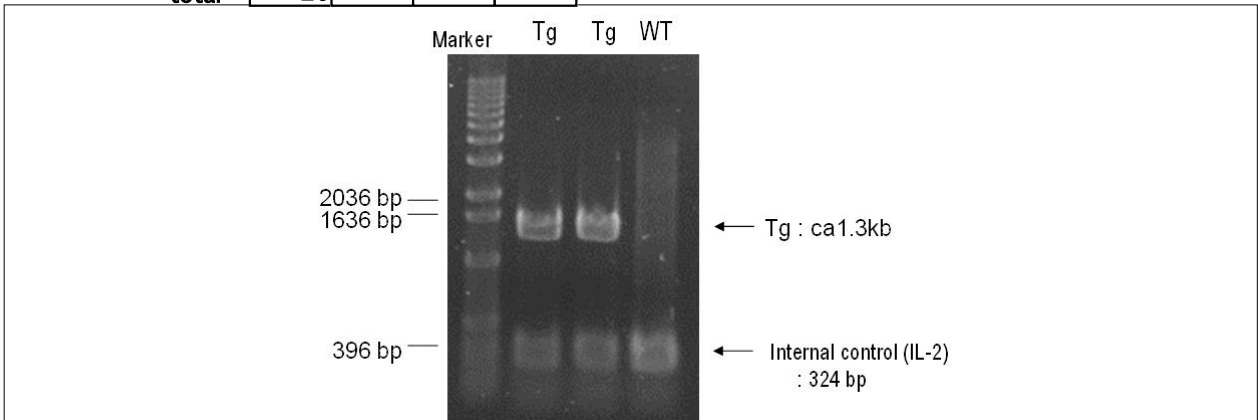

Marker
 Gene Ladder 100 (Wako, Cat. 316-06951)

Gel
 1% agarose in 1 X TAE (100 V, 30 min)

PCR products

primer_set	product_size	product
primer 1 primer 2	ca 1.3 kb	Tg
primer 3 primer 4	324 bp	Internal control (IL-2)

Comment