

RBRC 05654

Strain name B6.Cg-Tg(CAG-mRFP1,-SOX9,EGFP)1Haak

Primer	primer_name	primer_seq	length	
1	5CAG	CCTACAGCTCCTGGGCAACGTGC	23	RIKEN BRC
2		GTCGAGTGAGCTGTGTGTAGACG	23	RIKEN BRC
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol

Reaction Components


Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	68	300*	*Repeat for 30 cycles

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		0.6			
PCR buffer	10 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/μl	1.0			
primer 2	10 pmol/μl	1.0			
primer 3	10 pmol/μl	1.0			
primer 4	10 pmol/μl	1.0			
DNA		1.0			
Taq	1 U/μl	0.4			
total		20			

Standard PCR

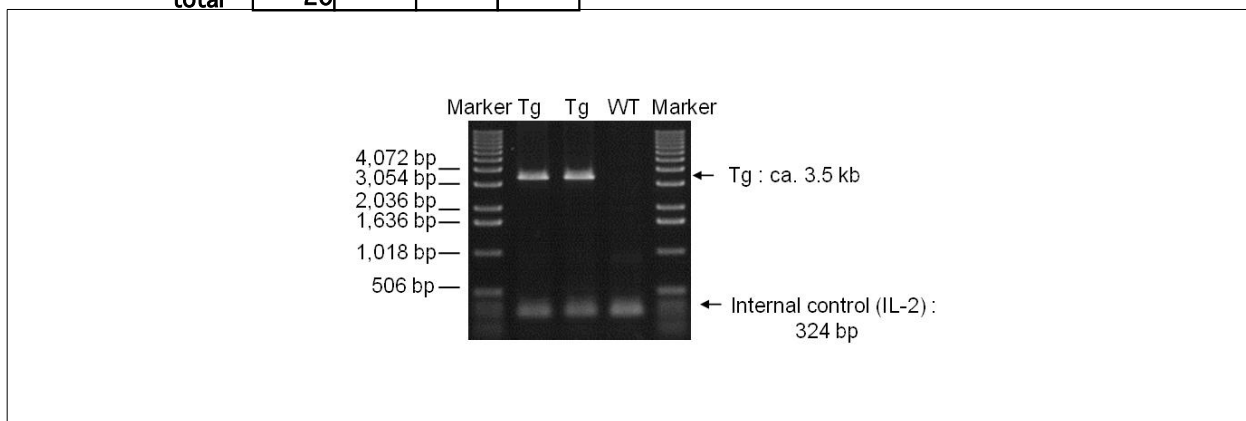

Marker
 1 kb DNA ladder (Invitrogen, Cat. 15615-024)
Gel
 1% agarose in 1 X TAE (100 V, 30 min)

PCR products			
primer_set	product_size	product	
primer 1	primer 2	ca. 3.5 kb	Tg
primer 3	primer 4	324 bp	Internal control (IL-2)

Comment Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.