

RBRC 05913

Strain name C57BL/6N-Tg(CAG-anti-Adipo)1Shb

Primer	primer_name	primer_seq	length	
1	5CAG	CCTACAGCTCCTGGGCAACGTGC	23	RIKEN BRC
2	mGBP191	ATCCTGGCCACAATGGCACA	20	Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol


### Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

### Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ $\mu$ l	1.0			
primer 2	10 pmol/ $\mu$ l	1.0			
primer 3	10 pmol/ $\mu$ l	1.0			
primer 4	10 pmol/ $\mu$ l	1.0			
DNA		1.0			
Taq	1 U/ $\mu$ l	0.2			
<b>total</b>		<b>20</b>			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

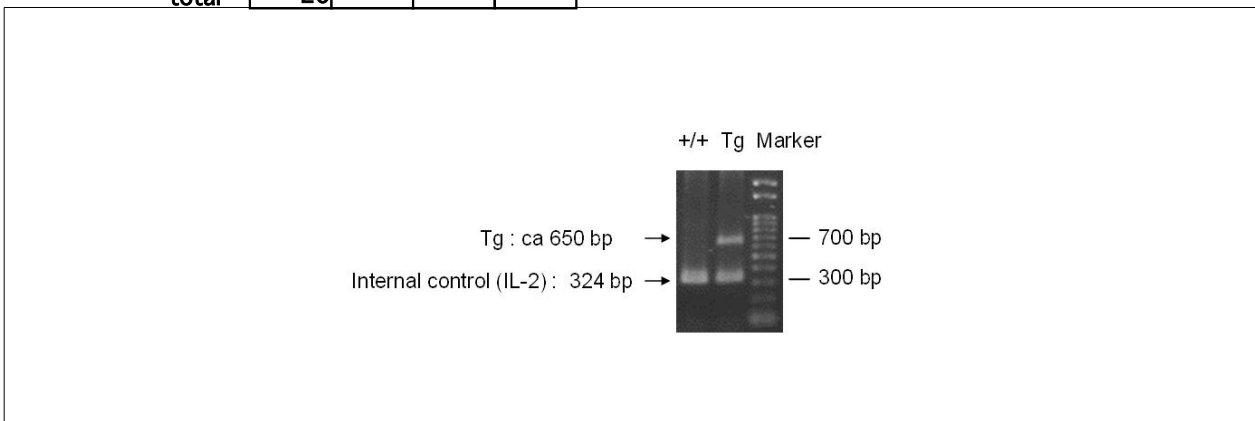

### Marker

Gene Ladder 100 (Wako, Cat. 316-06951)

### Gel

2% agarose in 1 X TAE (100 V, 30 min)

### PCR products

primer_set	product_size	product
primer 1 primer 2	ca 650 bp	Tg
primer 3 primer 4	324 bp	Internal control (IL-2)

**Comment** Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.