

RBRC 06569

Strain name B6(C3)-Ift80<tm1c(KOMP)Wtsi>

Primer	primer_name	primer_seq	length	Reference
1	IFT80-F6	GCTTTTGGTAAGGTGGGTCACACT	24	Depositor's primer
2	IFT80-R4	CACTCTCTGTCAGTGGTCCTTCA	24	Depositor's primer

Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)


Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Reagent	Conc	Reaction 1	2	3	4	(µL)
H2O		2.8				
PCR buffer	2 X	10.0				
dNTP	2.0 mM	4.0				
primer 1	10 pmol/µl	1.0				
primer 2	10 pmol/µl	1.0				
DNA		1.0				
Taq	1 U/µl	0.2				
total		20				

Marker Gene Ladder 100 (Wako, Cat. 316-06951)

Gel 2% agarose in 1 X TAE (100 V, 30 min)

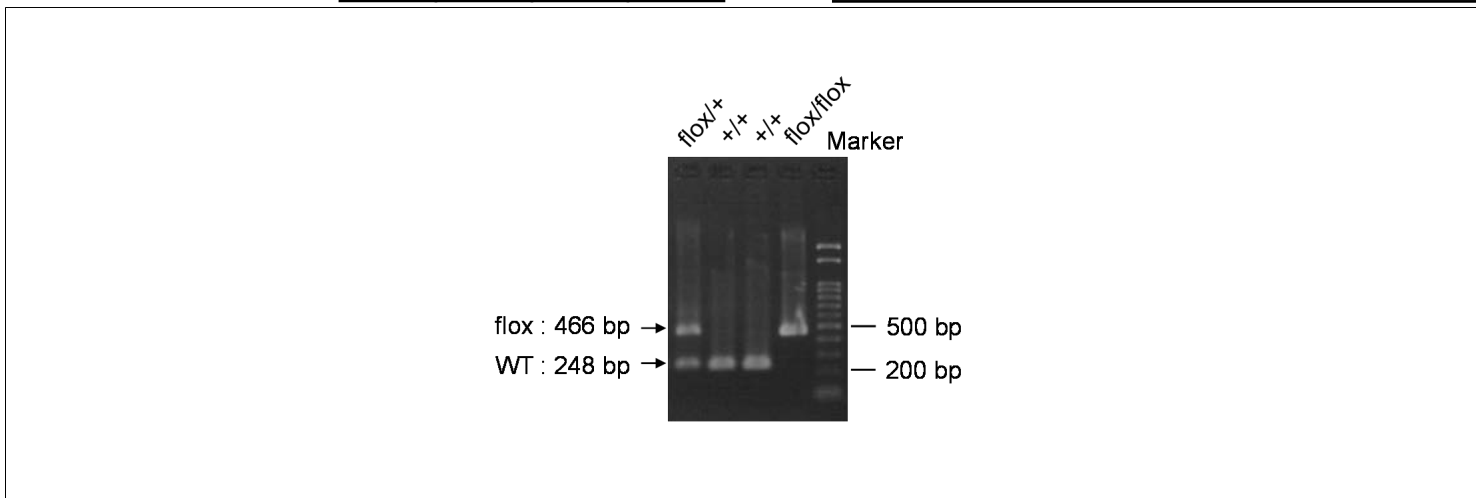

primer_set	product_size	product (allele type)	
primer 1	primer 2	248 bp	WT
primer 1	primer 2	466 bp	flox

Comment Standard PCR

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.