

RBRC 09303

Strain name C.129(B6)-Chst15<tm1Oh> (BALB/cCrSlc)

Primer	primer_name	primer_seq	length	Reference
1		GGCTATGGTGGTTCAAGCTGC	21	RIKEN BRC
2	Exon2R	GGCTCCTTGATGGCTGAGAATTTTCAC	26	Depositor's primer
3	NeoP	CTTCTATCGCCTTCTTGACGAGTT	24	Depositor's primer

Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)


Reagent	Conc	Reaction 1	2	3	4
H2O		1.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
primer 3	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

(μ L)

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Marker Gene Ladder 100 (Wako, Cat. 316-06951)
Gel 2% agarose in 1 X TAE (100 V, 30 min)

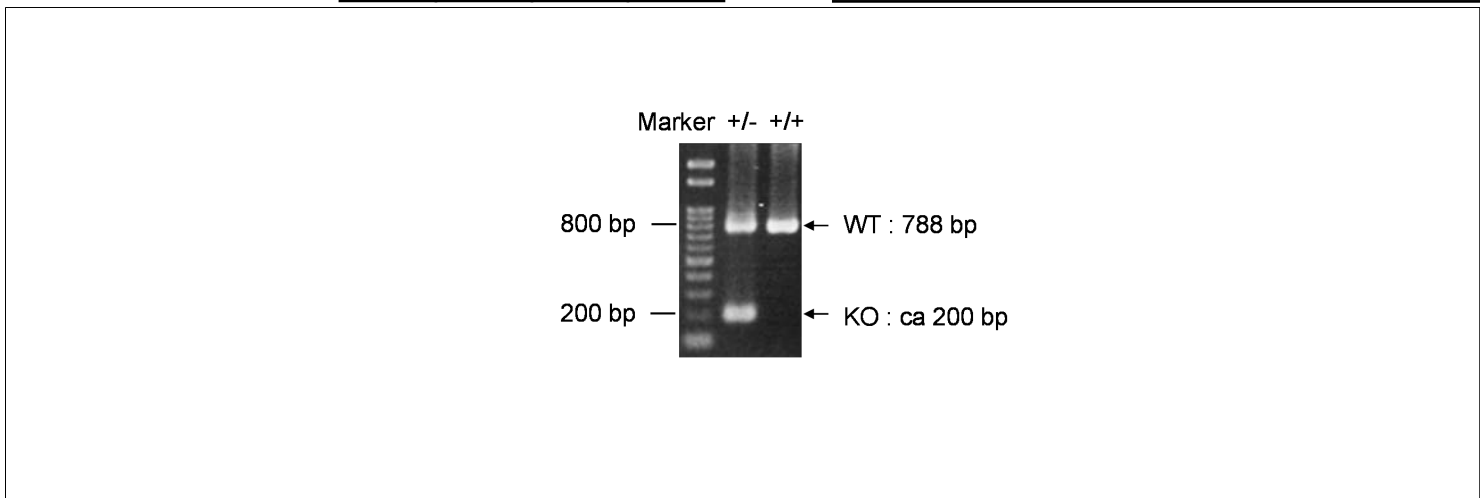

primer_set	product_size	product (allele type)
primer 1 / primer 2	788 bp	WT
primer 1 / primer 3	ca 200 bp	KO

Comment Standard PCR

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.