

RBRC 09415

Strain name ICR.Cg-Tg(Kiss1-AcGFP)3Nurep

Primer	primer_name	primer_seq	length	Reference
1	#72	CCTGGGTATCTCTCAAGTGCAGAAA	25	Depositor's primer
2	oIMR0872	AAGTTCATCTGCACCACCG	19	JAX Genotyping Protocol
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol

Reaction Components


Enzyme	TOYOBO KOD FX (Code.No : KFX-101)
--------	-----------------------------------

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Reagent	Conc	Reaction 1	2	3	4	(μ L)
H2O		0.8				
PCR buffer	2 X	10.0				
dNTP	2 mM	4.0				
primer 1	10 pmol/ μ l	1.0				
primer 2	10 pmol/ μ l	1.0				
primer 3	10 pmol/ μ l	1.0				
primer 4	10 pmol/ μ l	1.0				
DNA		1.0				
Taq	1 U/ μ l	0.2				
total		20				

Marker	Gene Ladder 100 (Wako, Cat. 316-06951)
Gel	2% agarose in 1 X TAE (100 V, 30 min)

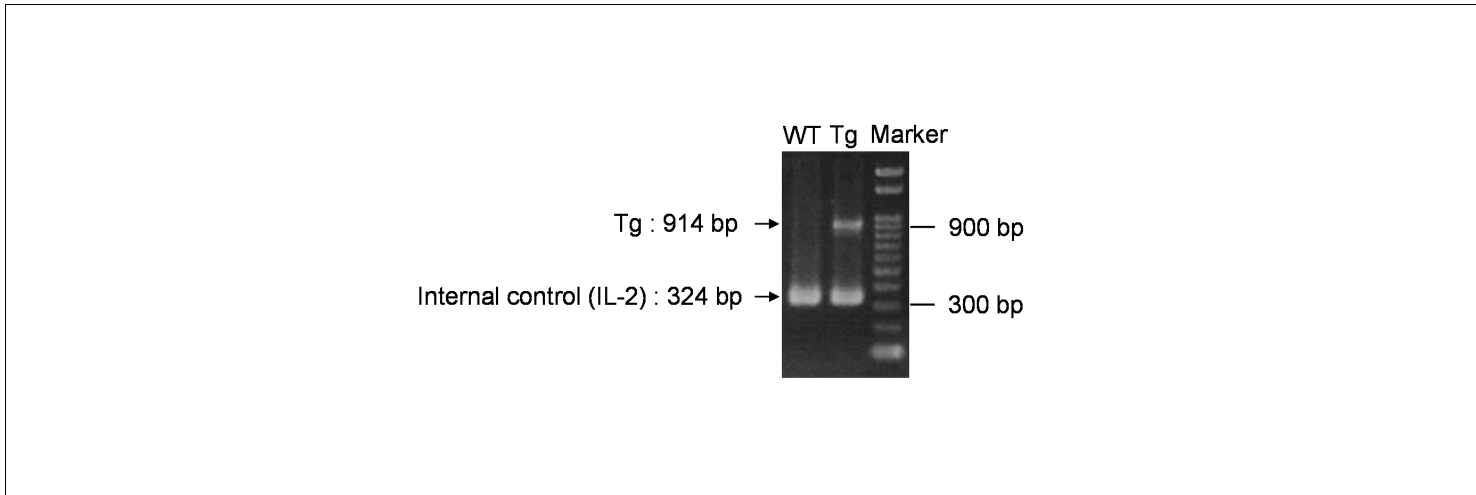

primer_set	product_size	product (allele type)
primer 1	914 bp	Tg
primer 3	324 bp	Internal control (IL-2)

Comment

--