

RBRC 10015

Strain name B6.Cg-Bcl6<tm1.1Toka>

Primer	primer_name	primer_seq	length	
1	Mh73_Bcl6YFP3	ACCAGCACACTCCAAGTAACA	21	Immunity, 34(6), 961-972 (2011).
2	Mh70_Bcl6YFPR	TGAACAGACGTTAAACTCCTATTA	25	Immunity, 34(6), 961-972 (2011).

Reaction Components


Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

	Conc	Run 1	Run 2	Run 3	Run 4
H2O		2.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/ μ l	1.0			
primer 2	10 pmol/ μ l	1.0			
DNA		1.0			
Taq	1 U/ μ l	0.2			
total		20			

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Standard PCR

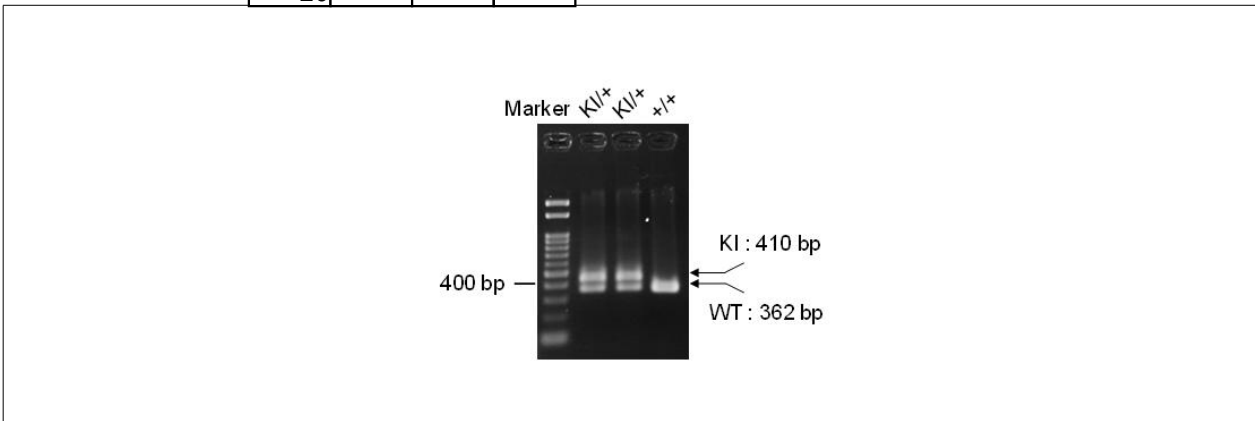

Marker
Gene Ladder 100 (Wako, Cat. 316-06951)

Gel
2% agarose in 1 X TAE (100 V, 30 min)

PCR products	primer_set	product_size	product
primer 1	primer 2	362 bp	WT
primer 1	primer 2	410 bp	KI

Comment

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.