

RBRC 10255

Strain name C.B6-Il17rb<tm1Tgi>

Primer	primer_name	primer_seq	length	Reference
1	il17rbWTf	TAACTGGCGTCCCGGGAGTG	20	Depositor's primer
2	il17rbWTr	GCTGGCAGGTTTTTCAGGCGC	20	Depositor's primer
3	neo_F	GCAGCGCATCGCCTTCTATCGC	22	RIKEN BRC

### Reaction Components


Enzyme TOYOBO KOD FX (Code.No : KFX-101)

### Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Reagent	Conc	Reaction 1	2	3	4	( $\mu$ L)
H2O		1.8				
PCR buffer	2 X	10.0				
dNTP	2.0 mM	4.0				
primer 1	10 pmol/ $\mu$ l	1.0				
primer 2	10 pmol/ $\mu$ l	1.0				
primer 3	10 pmol/ $\mu$ l	1.0				
DNA		1.0				
Taq	1 U/ $\mu$ l	0.2				
<b>total</b>		20				

**Marker** Gene Ladder 100 (Wako, Cat. 316-06951)  
**Gel** 2% agarose in 1 X TAE (100 V, 30 min)

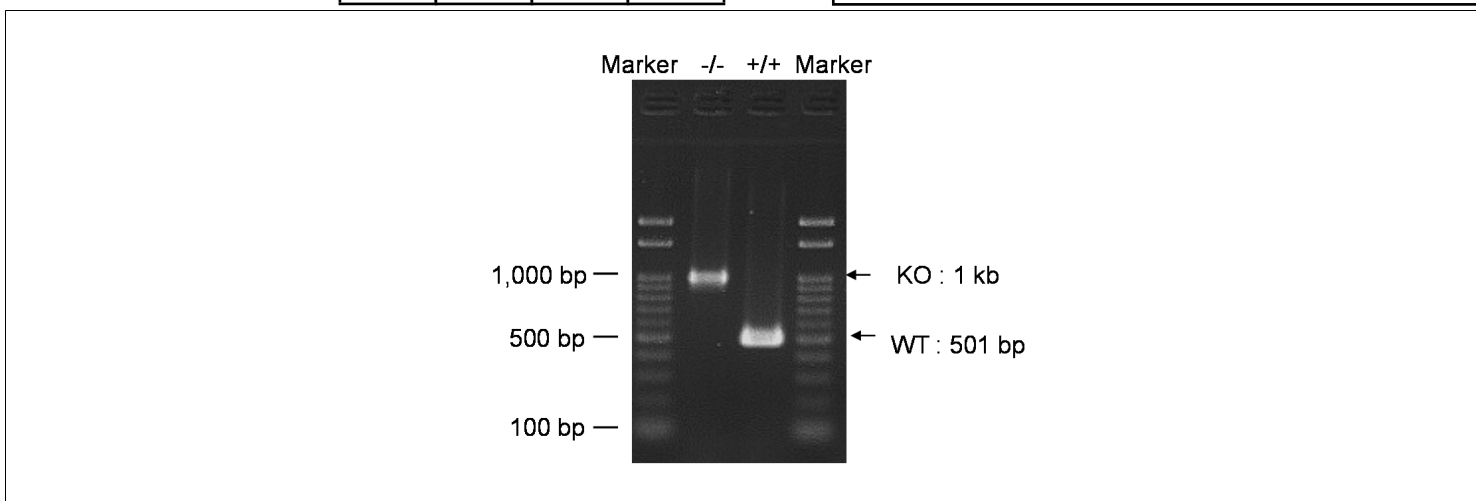

primer_set	product_size	product (allele type)	
primer 1	primer 2	501 bp	WT
primer 1	primer 3	1 kb	KO

+/+ (wild-type): 501 bp  
-/- (homozygous): 1 kb  
+/- (heterozygous): 501 bp and 1 kb

**Comment** Standard PCR

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.