

RBRC 10308

Strain name B6.Cg-Nfil3<tm1.2Kubo>

Primer	primer_name	primer_seq	length	Reference
1	E4bp4S1	TCAAAGTGGAGGCTTTGGAC	20	Depositor's primer
2	E4BP4KO AS1	CACAAGGACACCCAGACAGA	20	Depositor's primer

Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

step	temp (°C)	time (sec)	note
------	-----------	------------	------


Reagent	Conc	Reaction 1	2	3	4
H2O		2.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/μl	1.0			
primer 2	10 pmol/μl	1.0			
DNA		1.0			
Taq	1 U/μl	0.2			
total		20			

(μL)

1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Marker Gene Ladder 100 (Wako, Cat. 316-06951)

Gel 2% agarose in 1 X TAE (100 V, 30 min)

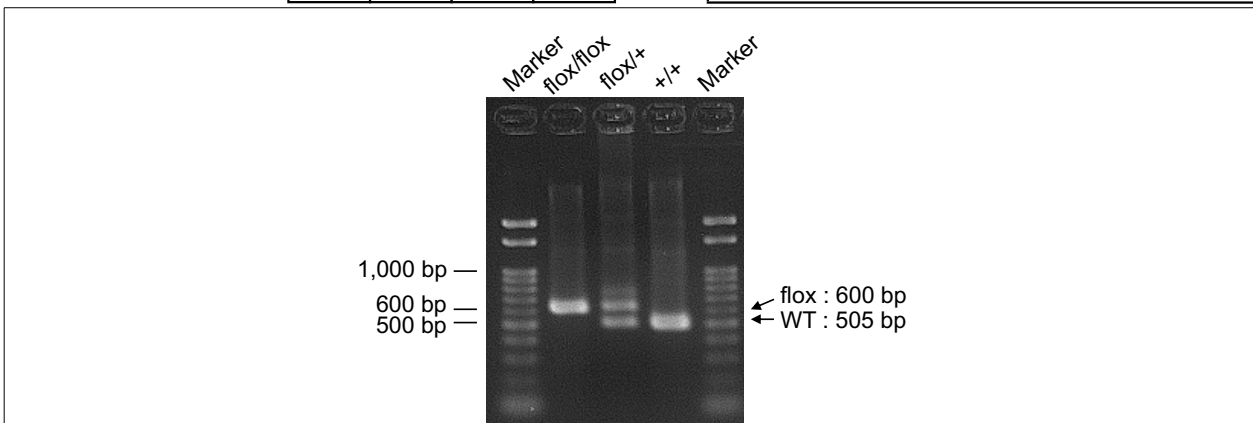

primer_set	product_size	product (allele type)
primer 1 primer 2	505 bp	WT
primer 1 primer 2	600 bp	flox

flox/flox : 600 bp
flox/+ : 600 bp, 505 bp
+/+ : 505 bp

Comment [Standard PCR](#)

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.