

RBRC 10387

Strain name C57BL/6-Pld2^{em1}Yaka

Primer	primer_name	primer_seq	length		Reference
1	PLD2-Cas9#1	GCCATCTATGACCTTCAGCCTCTGAAAGC	29		Depositor's primer
2	PLD2-Cas9#2	GTCTTCTGGCCCAACCCAGGCCTCCACTT	29		Depositor's primer

Reaction Components


Enzyme TOYOBO KOD FX (Code.No : KFX-101)

Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Reagent	Conc	Reaction 1	2	3	4	(μL)
H2O		2.8				
PCR buffer	2 X	10.0				
dNTP	2.0 mM	4.0				
primer 1	10 pmol/μl	1.0				
primer 2	10 pmol/μl	1.0				
DNA		1.0				
Taq	1 U/μl	0.2				
total		20				

Marker Gene Ladder 100 (Wako, Cat. 316-06951)
Gel 2% agarose in 1 X TAE (100 V, 30 min)

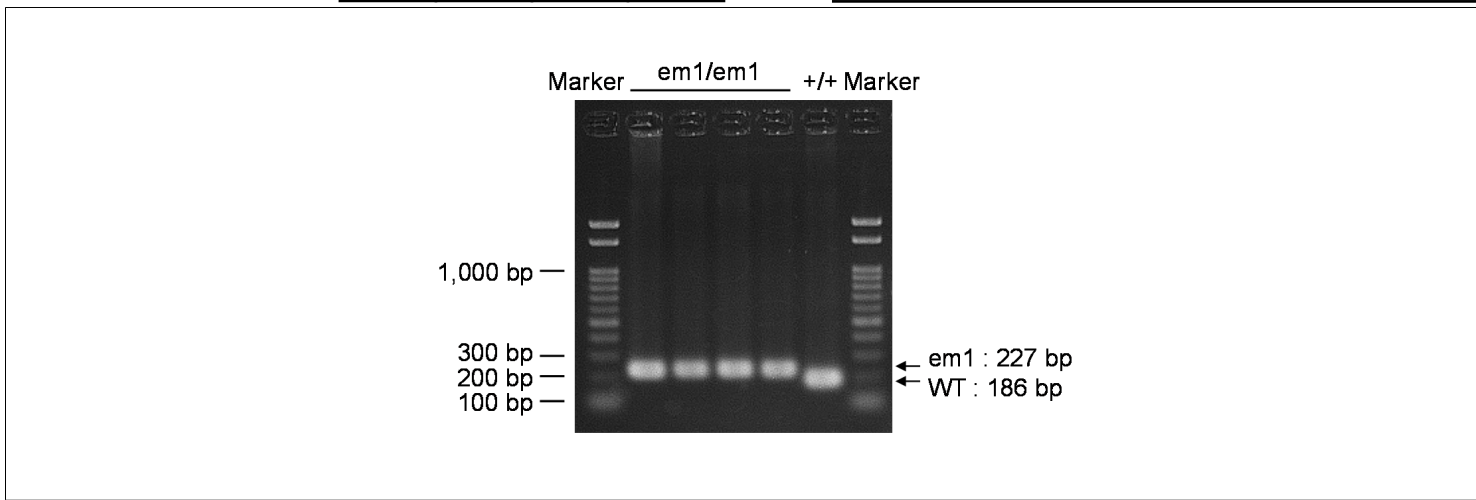

primer_set	product_size	product (allele type)	
primer 1	primer 2	186 bp	WT
primer 1	primer 2	227 bp	em1
			+/+ (wild-type): 186 bp em1/em1 (homozygous): 227 bp

Comment Standard PCR

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.