

RBRC 11845

Strain name B6.Cg-Tg(HCRT-EGFP)E2Stak

Primer	primer_name	primer_seq	length	Reference
1	hOX-promoter_F	GCAGCGGCCATTCCTTGG	18	Depositor's primer
2	EGFP32	TCAGCTCGATGCGGTTTAC	19	Depositor's primer
3	oIMR0042	CTAGGCCACAGAATTGAAAGATCT	24	JAX Genotyping Protocol
4	oIMR0043	GTAGGTGGAAATTCTAGCATCATCC	25	JAX Genotyping Protocol

Reaction Components

Enzyme TOYOBO KOD FX (Code.No : KFX-101)


Reaction Conditions

step	temp (°C)	time (sec)	note
1	94	120	
2	98	10*	
3	60	30*	
4	68	60*	*Repeat for 30 cycles

Reagent	Conc	Reaction 1	2	3	4
H2O		0.8			
PCR buffer	2 X	10.0			
dNTP	2.0 mM	4.0			
primer 1	10 pmol/μl	1.0			
primer 2	10 pmol/μl	1.0			
primer 3	10 pmol/μl	1.0			
primer 4	10 pmol/μl	1.0			
DNA		1.0			
Taq	1 U/μl	0.2			
total		20			

Marker Gene Ladder 100 (Wako, Cat. 316-06951)

Gel 2% agarose in 1 X TAE (100 V, 30 min)

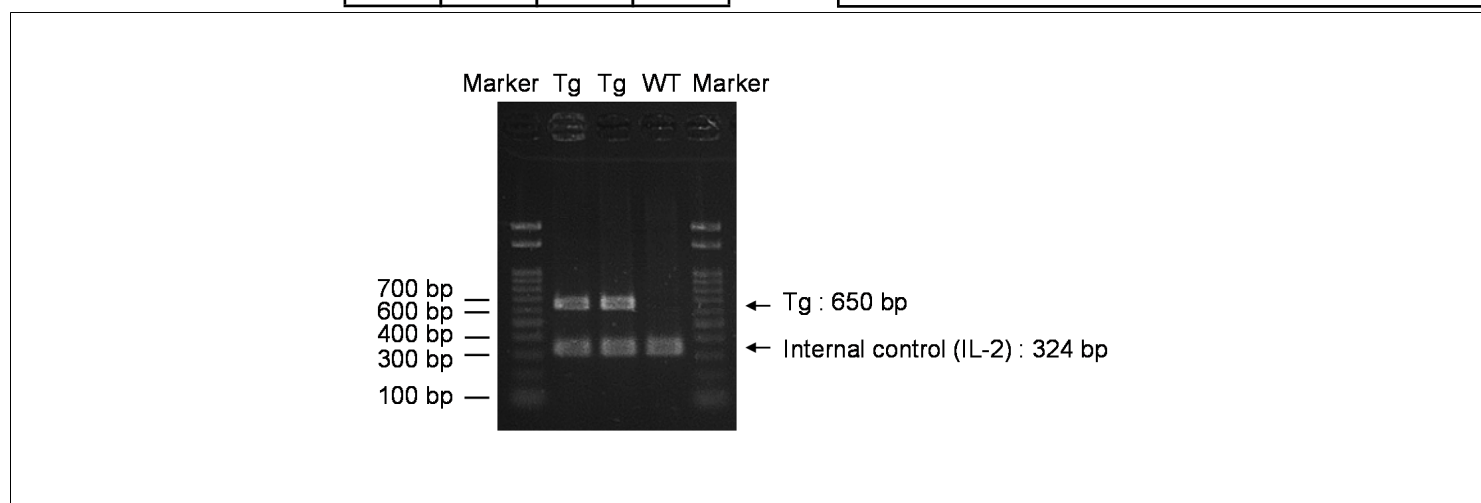

primer_set	product_size	product (allele type)
primer 1 primer 2	650 bp	Tg
primer 3 primer 4	324 bp	Internal control (IL-2)

Comment Standard PCR

Novagen KOD Xtreme™ Hot Start DNA Polymerase can be purchased from EMD.